500 kHz Lock-In Amplifier

SR860 — 500 kHz dual phase lock-in amplifier



- 1 mHz to 500 kHz frequency range
- Low-noise current and voltage inputs
- Touchscreen data display large numeric results, chart recordings, & FFT displays
- 10 MHz timebase input and output
- Dual reference mode
- GPIB, RS-232, Ethernet and USB
- HDMI video output

SR860 500 kHz Lock-In Amplifier

Superb performance. Outstanding value. They're what you've come to expect from a Stanford Research Systems lock-in amplifier. And they're delivered by the new SR860 500kHz Lock-in Amplifier, the latest in a line of innovative lock-ins from SRS. With unparalleled analog performance, sophisticated new digital signal processing features, a thoroughly modern, intuitive user interface, and a wide range of computer connectivity options, the SR860 is the ideal choice for any synchronous detection application.

With over 30 years of lock-in design experience SRS has made every effort to optimize each detail of the SR860. From a hefty toroidal transformer that eliminates switch-mode noise to iOS connectivity that brings the lock-in to your cell phone to advanced DSP filters that eliminate more noise while speeding up your experiment, the SR860 is truly the ultimate lock-in amplifier.

Signal Inputs

Lock-in performance starts at the front end. The SR860 front end offers both state-of-the art voltage and current input amplifiers. The voltage input is a switchable single-ended/ differential JFET-pair amplifier with 2.5 nV/ \sqrt{Hz} of noise at 1 kHz and under 10 nV/ \sqrt{Hz} of noise at 10 Hz. The voltage input has a 10 MΩ input impedance and can be AC or DC coupled. Input connector shields can be connected to the instrument ground through a user selectable 10 Ω (Ground) or 10 kΩ (Float) resistor.





SRS current source and several SRS preamplifiers

The SR860's built-in current amplifier represents a significant improvement over previous designs. The current input range is selectable from 1 μ A or 10 nA. The 1 μ A range has 400 kHz of bandwidth and 130 fA/ \sqrt{Hz} of noise, while the 10 nA range offers 2 kHz of bandwidth and 13 fA/ \sqrt{Hz} of noise.

While the built in voltage and current amplifiers are suitable for most applications, the SR860 is also compatible with the complete range of specialized pre-amplifiers offered by SRS. The SR550 (FET input), SR552 (BJT input), SR554 (transformer input), SR555 (120 kHz current amp), and SR556 (low-noise current amp) can all be powered directly from the SR860's rear-panel preamp power port. The SRS website shows you how each of these preamps can be used to optimize measurements in a wide variety of experimental situations.

Sensitivity and Input Range

As with previous instruments, the Sensitivity setting of the SR860 is the voltage (or current) which produces a full scale output. But unlike previous designs, the input range of the SR860 can be explicitly set from the front panel without having to consult a confusing "dynamic reserve" equation. Simply choose the sensitivity required by your experiment and then select the smallest input range that doesn't overload. That's it.

The SR860's effective dynamic reserve is simply the ratio of these two settings. For instance with a 10 nV sensitivity setting and a 300 mV input range the effective dynamic reserve of the SR860 is 3×10^7 , or nearly 150 dB.



Large numeric readout with strip chart recording



Output Time Constants and Filtering

The SR860 offers traditional RC-response output time constants from 1 μ s to 30 ks with rolloffs of 6, 12, 18, and 24 dB/oct. But in addition, the SR860 also offers advanced digital filters which can significantly reduce measurement time while increasing signal to noise. Below 3 s, the advanced filters are Gaussian FIR filters which at the same noise bandwidth as an RC filter have significantly better rise-time and stopband attenuation. These filters also have symmetric rise and fall profiles which preserves feature shapes while scanning in frequency. At time constants longer than 3 s, the advanced filters are linear phase IIR filters which settle nearly twice as fast as their RC counterparts for equivalent stopband attenuation.

Synchronous filtering may also be selected at reference frequencies below 4 kHz. The synchronous filter notches out multiples of the reference frequency and is extremely useful in making low frequency measurement where multiples of the reference frequency would otherwise show up in the lock-in output. Unlike previous designs the synchronous filter in the SR860 can be selected with no loss of output resolution.



SR540 Optical Chopper

Reference Channel


The SR860 has a specified reference frequency range of 1 mHz to 500 kHz. Detection can be done at the fundamental of the reference frequency, or at up to the 99th harmonic. Several reference modes are available: Internal mode uses the SR860's precision internal oscillator as the reference. External mode locks to an external sine or TTL signal. In Dual Mode, the lock-in detects at the difference frequency between the internally set reference frequency and an externally applied sync signal allowing direct recovery of a double-modulated signal. Finally, in Chop mode, the SR860 provides a digital PID (Proportional-Integral-Derivative) controller signal to synchronize an SR540 Optical Chopper to the internal oscillator of the lock-in. By having the lock-in directly control the chopper, frequency drift can be virtually eliminated.

Sine Output

The SR860 offers a precision sine wave output which can be set with 6 digits of frequency resolution and an amplitude range from 1 nV to 2 V. The SR860 output is unique in that it can be configured as a single-ended or as a differential (balanced) signal. A DC offset of up to \pm 5 V can be applied to the sine output. A rear-panel logic-level sync signal synchronized to the sine output is also provided.

Timebase

Rear-panel 10 MHz inputs and outputs are provided allowing the SR860 to be locked to an external frequency reference (such as the FS725 10 MHz Rubidium Frequency Standard). Alternatively, the 10 MHz output from the SR860 can be used to synchronize several lock-ins or other test equipment with a 10 MHz timebase input.

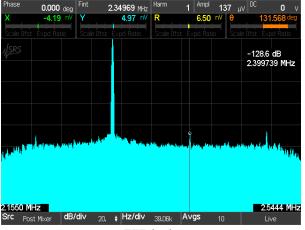

FFT display

FFT Displays

Lock-in amplifiers are traditionally time-domain instruments but sometimes it's easier to understand a signal in the frequency domain. The SR860 is at home in both worlds. An FFT display shows the spectrum of the input signal at the front end, the post-mixer signal, or the spectrum of the signal after the time-constant filters. Using the FFT displays can simplify tracking down sneaky noise sources that would otherwise get lost in the "one-number" output of a conventional lock-in amplifier.



Trend analysis

Front-Panel Touchscreen Display

The center of the SR860's front panel is a full-color $640 \times$ 480 touchscreen which can be set to display up to 4 channels of data. (When a dark lab is required, the LCD screen can be blanked from the front panel or from the remote interface.) Each data channel can be configured to display X, Y, R, Θ , Aux In (1-4), Aux Out (1-2), X noise, Y noise, Sine Amplitude, Sine Out DC Level, Reference Phase, or Reference Frequency. The screen can be set up to show the data channels as large numbers, easily visible from across the room, or as a "strip-chart" display showing a complete history of each channel with selectable time scales from 0.5 s/div to 2 days/div. Even when not displaying a measurement, all measurements are always being saved by the SR860 insuring that no data is ever lost. The touchscreen also continually displays key lockin setup parameters such as phase, reference frequency, and sine amplitude. And a rear-panel HDMI port allows the LCD screen to be viewed on any HDMI monitor or TV.

It should be emphasized that as useful as the touchscreen display is for displaying data, it's not necessary to use the touchscreen to control common instrument functions on the SR860. All commonly used controls: time constant, reference frequency, sine amplitude and offset, input configuration, and more, are controllable with dedicated front panel knobs or buttons. Infrequently accessed configuration settings, such as the TCP/IP settings and other communication settings, are accessed through menus shown on the touchscreen.

Computer Connectivity

The SR860 comes standard with virtually every remote interface imaginable. GPIB (IEEE488.2) and RS-232 are of course provided, as well as USB (Test and Measurement Class) and Ethernet (VXI-11 and telnet). The SR860 hosts its own webserver allowing the instrument to be monitored and controlled remotely with just a browser. And when connected to a Wi-Fi enabled network, the SR860 can be controlled remotely from an iPhone or iPad with free iOS app downloadable from the Apple App Store. The app allows remote configuration of the instrument and monitoring of live data from your handheld device.

A front-panel USB port allows data and screen-shots to a USB flash drive. Data can be saved either as comma delimited files or MATLAB compatible .MAT files. Incorporating screen shots and data into reports or spreadsheets has never been easier.





SR860 Specifications

Single-ended or differential

1 fA to $1 \mu A$ (current input)

1 V (peak) or 1 µA (peak)

 $10 M\Omega + 25 pF$, AC or DC coupled

Signal amplitude must be less than

-80 dB (<100 kHz), -60 dB (>100 kHz)

 $1 k\Omega$ or 100Ω to virtual ground

10 mV to 1 V (peak)

1 µA or 10 nA (peak)

 $\pm 1\%$ (<200 kHz)

120 dB (typ)

30% of input range 2.5 nV/√Hz at 1 kHz

0.001 Hz to 500 kHz

<0.001° rms at 1 kHz (100 ms, 12 dB/oct)

<0.02°/°C below 200 kHz

<0.2°/°C below 2 MHz

(10 mV input range, typ.)

90 dB to 1 kHz (DC coupled),

10 MHz In/Out (phase locks the internal frequency to other SR860s)

<0.002°/°C below 20 kHz (DC coupled)

Detect at N× f_{ref} (N<99 and (N× f_{ref})<2 MHz) Detect at $f_{dual} = |f_{int} - f_{ext}|$ SR860 drives SR540 Chopper (via

Aux Out 4) to lock the chopper to fint

Signal Channel

Voltage inputs

Sensitivity (output scale) 1 nV to 1 V (voltage input) Voltage input range Current input range Max input Input impedance Voltage input Current input Gain accuracy Noise (rms)

Harmonic distortion CMRR Dynamic reserve

Reference Channel

Frequency range Timebase

Input impedance $1 M\Omega$ or 50Ω Phase setting resolution 360/2³² degrees Phase noise <0.0001° rms at 1 kHz Int. ref (100 ms, 12 dB/oct)

Ext. ref (typ)

Phase drift (typ)

Harmonic detection

Dual F reference Chopper reference

Demodulator

DC stability Time constants Low pass filters

Filter slope Synchronous filter Harmonic rejection Low latency output

Digital output values have no drift 1 µs to 30ks Typical RC-type filters or advanced Gaussian/Phase-Linear filters 6, 12, 18 or 24 dB/oct rolloff Available below 4 kHz -80 dB Rear-panel BlazeX output with <2 µs delay (plus LPF rise/fall times).

Internal Oscillator

0.001 Hz to 500 kHz Frequency range Frequency accuracy $25 \text{ ppm} + 30 \,\mu\text{Hz}$ (with internal timebase) External timebase 10 MHz timebase input/output (on rear panel) Frequency resolution 6 digits or 0.1 mHz, whichever is greater

Sine Output

Sille Output	
Outputs Output impedance Amplitude Amplitude resolution	Single-ended or differential 50Ω source 1 nVrms to 2 Vrms (amplitude is differential into 50Ω loads). Output is halved when using in single-ended mode. Output is doubled when driving a high impedance load. 3 digits or 1 nV, whichever is greater
DC offset Offset resolution Output limit	± 5 V, differential or common mode 3 digits or 0.1 mV, whichever is greater ± 6 V, sum of DC offset and peak amplitude
Sync	Logic level sync on rear panel (via BlazeX output)
Data	
Data channels	Four data channels are displayed and
Data sources	graphed (green, blue, yellow, orange) Each data channel can be assigned any of these data sources: X, Y, R, θ , Aux In 1 to 4, Aux Out 1 to 2, X noise, Y noise, Sine Out Amplitude, Sine Out DC Level, Reference Phase, or
Data history	Reference Frequency. All data sources are continually stored at all chart display time scales. The complete stored history of any data source can be displayed sat any time.
Offset	X, Y and R may be offset up to $\pm 999\%$
Expand	of the output scale X, Y and R may be expanded by $\times 10$ or $\times 100$
Ratio	X and Y may be ratioed by Aux In 3. R may be ratioed by Aux In 4.
Capture buffer	I Mpoint internal data storage. Store (X), (X and Y), (R and θ), or (X, Y, R and θ) at sample rates up to 1.25 MHz. This is in addition to the data histories for the chart display.
Data streaming	Realtime streaming of data, either (X), (X and Y), (R and θ), or (X, Y, R and θ) at sample rates up to 1.25 MHz over the Ethernet interface.
Scanning	One of the following parameters may be scanned: f_{int} , Sine Out Amplitude, Sine Out DC Level, Aux Out 1 or 2
FFT	

FF

Source

Record length Averaging

Input ADC, demodulator output, or filter output 1024 bins Exponential rms



Inputs and Outputs	
CH1 output	Proportional to X or R $(\pm 10 \text{ V full scale thru } 50 \Omega)$
CH2 output	Proportional to Y and θ (±10 V full scale thru 50 Ω)
X and Y outputs	Proportional to X and Y (rear panel) $(\pm 10 \text{ V full scale thru } 50 \Omega)$
BlazeX	Low latency ouput of X, ± 2.0 V full scale or logic level reference sync
Aux outputs	output, either thru 50Ω 4 BNC D/A outputs, ± 10.5 V thru 50Ω ,
Aux inputs	1 mV resolution 4 BNC A/D inputs, ±10.5 V,
Trigger input	1 mV resolution, $1M\Omega$ input TTL input triggers storage into the internal capture buffer
Signal monitor	Analog output of the signal amplifier
HDMI	Video output to external monitor or TV (640 × 480, 60 Hz)
Timebase I/O	1 Vrms, 10 MHz clock to synchronize internal reference to other units

General

Interfaces	GPIB (IEEE-488.2), RS-232, USB
	and Ethernet
USB flash	Front-panel slot for USB flash storage
	of screen shots and data
Preamp power	9-pin D connector to power
1 1	SRS preamps
Power	60 W, 100/120/220/240 VAC, 50/60 Hz
Dimensions	17"×5.25"×17" (WHD)
Weight	22 lbs.
Warranty	One year parts and labor on defects
	in materials and workmanship

Ordering Information

SR860	500 kHz lock-in amplifier
SR550	Voltage preamplifier
	$(100 \mathrm{M}\Omega, 3.6 \mathrm{nV}/\sqrt{\mathrm{Hz}})$
SR552	Voltage preamplifier
	$(100 \mathrm{k}\Omega, 1.4 \mathrm{nV}/\sqrt{\mathrm{Hz}})$
SR554	Transformer preamplifier
	$(0.091 \text{nV}/\sqrt{\text{Hz}})$
SR555	Current preamplifier
SR556	Current preamplifier
SR540	Optical chopper



SR860 front panel



SR860 rear panel